

查询"L10671W"供应商

NEW

COMPACT UV-VIS S2D2 FIBER LIGHT SOURCE L10671

OVERVIEW

The L10671 is the compact UV-VIS fiber light source using a S2D2 (Stable and Small Deuterium) lamp. The L10671 is the light source that offers high radiant power and high stability despite its compact size. Besides being compact and easy to carry, the L10671 operates on a low voltage. These features make the L10671 ideal for all types of portable devices requiring a UV-VIS light source.



* Light guide is sold separately.

FEATURES


- Compact Size: 72 mm × 40 mm × 90 mm
- High Stability: Fluctuation 0.004 % p-p (Typ.)
(Equivalent to 2×10^{-5} A.U.)
- External Control
- Shutter Function

APPLICATIONS


- Spectrophotometry
- Environmental Measurement
- Absorption Spectrum Measurement
- Micro Titer Plate Reader
- Pharmaceutical Testing
- High-Performance Liquid Chromatography
- Film Thickness Measurement
- Biological Measurement
- Semiconductor Inspection
- Optical Component Inspection

LIGHT OUTPUT STABILITY

● at 230 nm

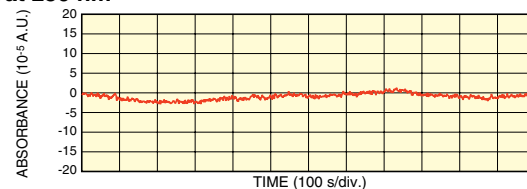


● at 600 nm

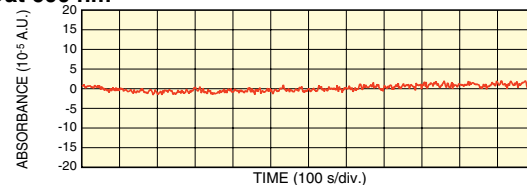

TLSOB0107EA

SPECTRAL DISTRIBUTION (TYPICAL DATA)

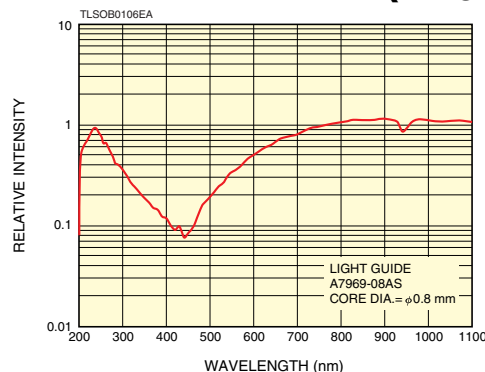

HAMAMATSU

SPECIFICATIONS

查询"L10671W"供应商

Parameter		Specifications
Spectral Distribution		200 nm to 1600 nm
Built-in Lamp		S2D2 Lamp 0.5 W Tungsten Lamp
Output Stability	Fluctuation (Typ.)	0.004 % p-p (Equivalent to 2×10^{-5} A.U.)
	Drift (Max.)	± 0.25 %/h
Input Voltage (AC) ①②		100 V to 240 V (50 Hz/60 Hz)
Power Consumption ①		10 VA
Weight		200 g
Light Guide Connector		SMA905D
Shutter		Included
External Control	Lamp ON/OFF	
	Shutter open/close	
	Status Signal for Lamp ON/OFF and Main Power ON/OFF	
Main Unit Switch		Shutter open/close
Operating Ambient Temperature		0 °C to +40 °C
Operating Ambient Humidity		Below 80 %
Storage Temperature		-10 °C to +60 °C
Storage Humidity		Below 85 %
Applicable Standard		CE Marking

NOTE: ①When used with AC adapter (supplied).

②Input voltage to light source unit is 12 V dc (8.5 V dc to 13.2 V dc).

REPLACEMENT LAMP

Type No.	Lamp	Lamp Life ③
L10671D	S2D2 Lamp	800 h
L10671W	Tungsten Lamp	1000 h

NOTE: ③Guaranteed life. Defined as the time when the radiant intensity falls to 50% of its initial value or the radiant output fluctuation exceeds 0.05 %.


LIGHT GUIDE ④ (SOLD SEPARATELY)

Type No.	Core Diameter	Length
A7969-06AS ⑤	0.6 mm	1.0 m
A7969-08AS ⑤	0.8 mm	


NOTE: ④Other shapes and combinations are available upon request.

⑤AS type is a UV-resistant light guide.

UV RESISTANT CHARACTERISTICS (at 215 nm)

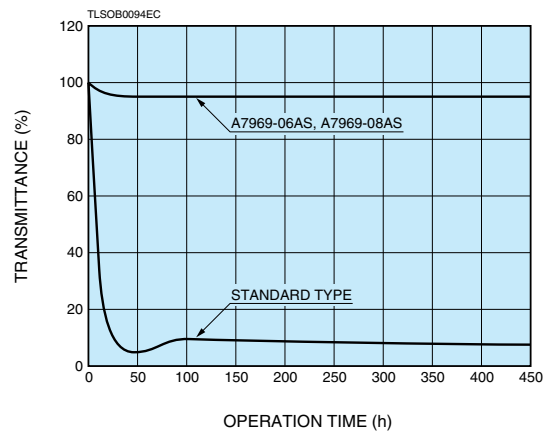

DIMENSIONAL OUTLINES (Unit: mm)

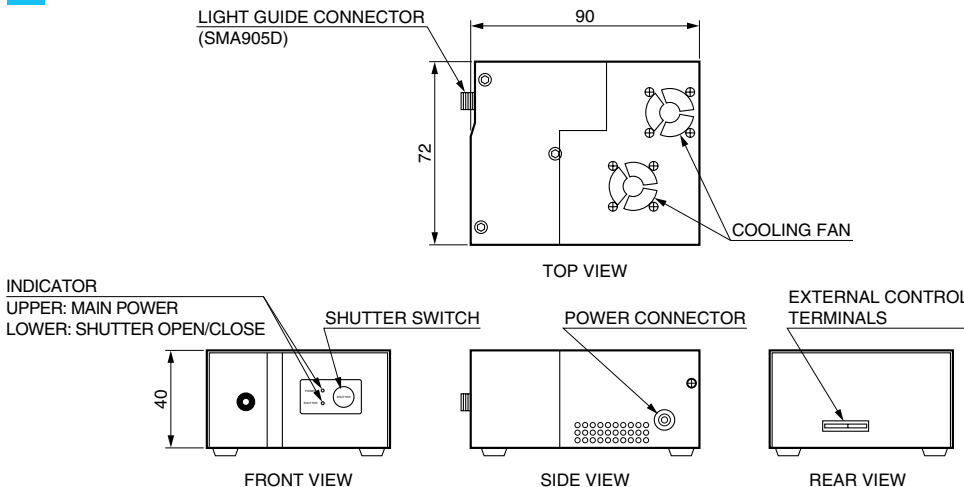

ACCESSORIES (Supplied)

- AC Adapter
- External Control Cable
- Instruction Manual

EXTERNAL CONTROL TERMINALS PINOUT

PIN No.	Signal Name
①	Deuterium lamp ON/OFF
②	Tungsten lamp ON/OFF
③	Shutter open/close
④	Deuterium lamp status signal
⑤	Tungsten lamp status signal
⑥	No connection
⑦	Main power ON/OFF status signal
⑧	No connection
⑨	No connection
⑩	GND

TLSOA0104EA

CE This product complies with safety standard IEC61010-1, EMC standard EN61326.

Subject to local technical requirements and regulations, availability of products included in this promotional material may vary. Please consult with our sales office. Information furnished by HAMAMATSU is believed to be reliable. However, no responsibility is assumed for possible inaccuracies or omissions. Specifications are subject to change without notice. No patent rights are granted to any of the circuits described herein. ©2007 Hamamatsu Photonics K.K.

HAMAMATSU

WEB SITE www.hamamatsu.com

HAMAMATSU PHOTONICS K.K., Electron Tube Division

314-5, Shimokanzo, Iwata City, Shizuoka Pref., 438-0193, Japan, Telephone: (81)539/62-5248, Fax: (81)539/62-2205

U.S.A.: Hamamatsu Corporation: 360 Foothill Road, P. O. Box 6910, Bridgewater, N.J. 08807-0910, U.S.A., Telephone: (1)908-231-0960, Fax: (1)908-231-1218 E-mail: usa@hamamatsu.com

Germany: Hamamatsu Photonics Deutschland GmbH: Arzbergerstr. 10, D-82211 Herrsching am Ammersee, Germany, Telephone: (49)8152-375-0, Fax: (49)8152-2658 E-mail: info@hamamatsu.de

France: Hamamatsu Photonics France S.A.R.L.: 19, Rue du Saule Trapu, Parc du Moulin de Massy, 91882 Massy Cedex, France, Telephone: (33)1 69 53 71 00, Fax: (33)1 69 53 71 10 E-mail: infos@hamamatsu.fr

United Kingdom: Hamamatsu Photonics UK Limited: 2 Howard Court, 10 Tewin Road Welwyn Garden City Hertfordshire AL7 1BW, United Kingdom, Telephone: 44-(0)1707-294888, Fax: 44-(0)1707-325777 E-mail: info@hamamatsu.co.uk

North Europe: Hamamatsu Photonics Norden AB: Smidesvägen 12, SE-171-41 SOLNA, Sweden, Telephone: (46)8-509-031-00, Fax: (46)8-509-031-01 E-mail: info@hamamatsu.se

Italy: Hamamatsu Photonics Italia: S.R.L.: Strada della Moia, 1/E, 20020 Arese, (Milano), Italy, Telephone: (39)02-935 81 733, Fax: (39)02-935 81 741 E-mail: info@hamamatsu.it

TLSO1067E01
SEPT. 2007 IP